

NEUATION



A WIN-WIN PROPOSITION

iTherm CH25
(Cooling - Heating Incubator)

iTherm CH25

(Cooling - Heating Incubator)

Cell and tissue culture work requires a safe, reliable and contaminant free surroundings temperature control from 15°C below ambient to 60°C, & the lowest achievable temperature is up to 2°C. Yes, there are many alternates available but when there is a lot at stake, you cannot rely on shoddy equipment that do not guarantee consistent and reliable temperature for however long you want. Therein the immaculate quality and design of Neuation will help you. With Neuation's iTherm CH25 you get a steady and reliable temperature with a wide range of control from 15°C below ambient to 60°C. It comes with an internal power outlet for a small shaker or mixer and can hold various volumes up to 2L. This incubator is designed to deliver quality results with pinpoint precision.



Versatile Perfection



Temperature Control Range
from 15°C below ambient to 60°C



Temp. Accuracy of
0.50C (at 370C)



Large Capacity Of
22 Litres, Allows To
Store Large Batches



Internal Power Outlet
For Installing Other
Instruments



Comes with an inbuilt
calibration function
of up to $\pm 5^{\circ}\text{C}$

Save Sample From
Cross Contamination



Easy To Move And Use At
Different Locations




Rapid Heating And Uniform
Temperature Across The Block
Throughout Its Process



Brings Down The
Temperature And Eliminates
The Need For A Compressor



iTherm CH25 Temperature Range Visualization

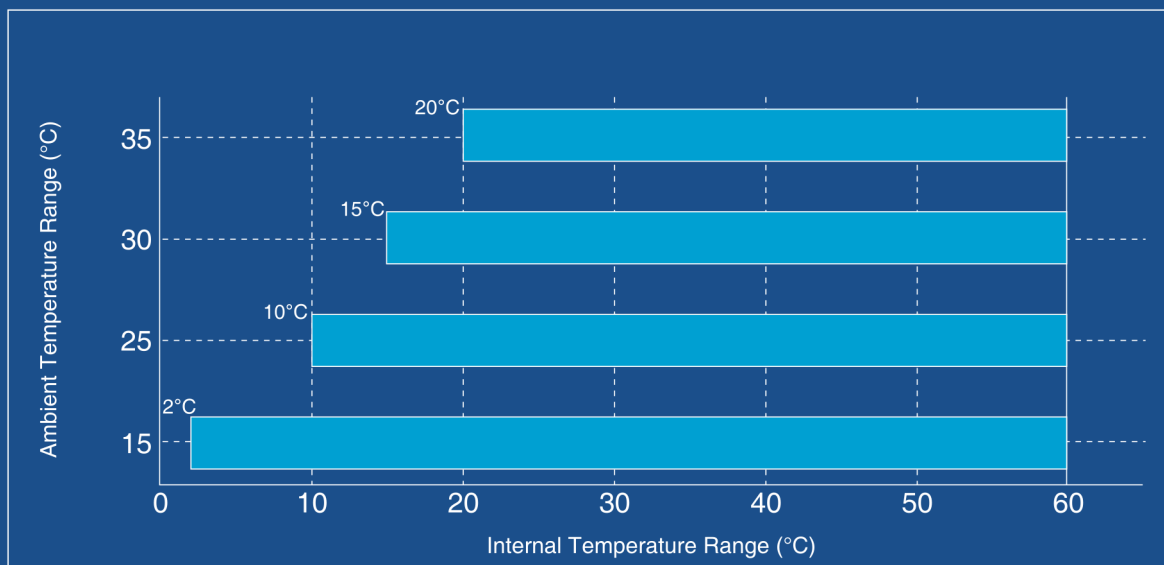

The iTherm CH25 provides temperature control from 15°C below ambient to 60°C, & the lowest achievable temperature is up to 2°C.

- Ambient 25°C » Range = 10°C to 60°C
- Ambient 15°C » Range = 2°C to 60°C

Specifications

iTherm CH25	
Specifications	Details
Temperature Range	15°C below Ambient to 60°C
Lowest Achievable Temperature	2°C
Temperature Accuracy	0.5°C (at 37°C)
Temperature Uniformity	+/- 1.5°C (at 37°C)
Temperature Increments	+/- 1°C
Exterior Dimension	410 x 340 x 480 mm
Interior Dimension	260 x 235 x 360 mm
Capacity	22L
Weight	7.4 Kgs
Electrical	100-240V, 50-60Hz

Ordering information

iTherm CH25		
Product Name	Cat #	
iTherm CH25	IN01-2232 EU/UK/AUS/US	Incubator With Ambient -15°C to 60°C

Caresource For Scientific & Laboratory Equipment Trading LLC
www.caresourceglobal.com

AN ISO 9001 : 2015 & ISO 13485 : 2016 CERTIFIED MFG. COMPANY

*We also execute OEM projects. Specifications can be changed without notice for quality improvement.